


Impeccable
Credentials



CONTECH[®]
Instruments Ltd.

Where innovation never stops...

DENSITY DETERMINATION OF PLASTIC GRANULES by CONTECH[®] Density Balance

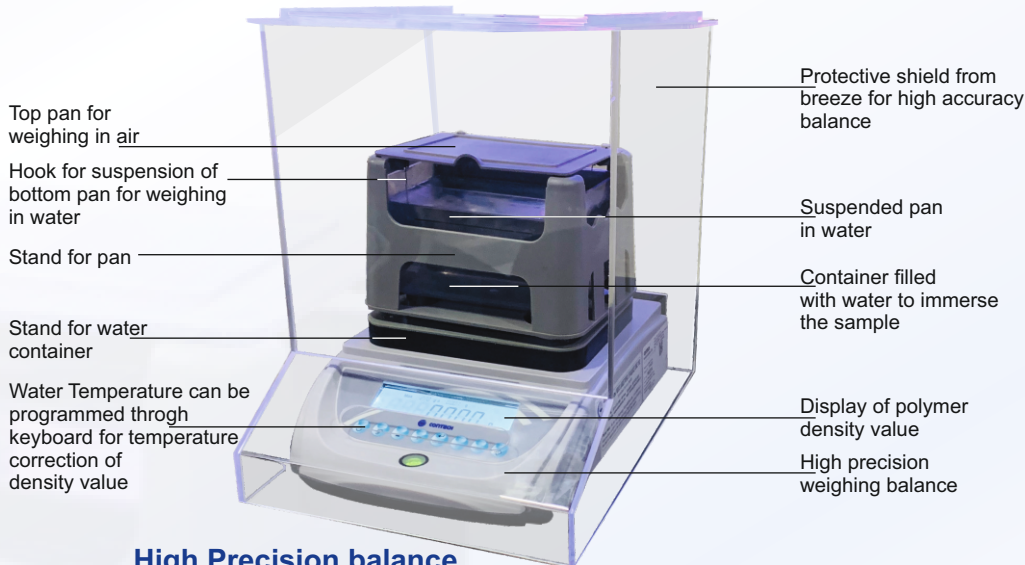

High Precision balance accuracy 1mg

Density Balance for Plastic Industries based on Archimedes principle.

- Density determination kit for plastic granules.

- > Granules in container weighed in air.
- > Granules transferred to a stainless steel fine mesh container with lid & weighed in water.
- > The software calculates polymer density & displays directly on the balance.
- > Results can be printed in the printer attached to the balance.



Stainless steel
fine mesh spherical
sample container



Water Container



Printer
Print out of density on the printer
attached to the balance



CONTECH INSTRUMENTS LTD.

301, Punit Indl. Premises, Turbhe Naka, Navi Mumbai - 400 705.

Tel. : +91 22 6139 3000 (100 lines)

E-mail : sales@contechindia.in

Where innovation never stops.....
www.contechindia.com