


FLAMEPROOF BALANCES


0.1 mg - 5 kg

CONTECH High Precision Balances incorporate advanced Electro Magnetic Force Compensation (EMFC) technology for offering unmatched accuracy and reliability, coupled with unique ADC circuit with inbuilt temperature compensation to offer a very high stability upto a resolution of 6 million counts. All adjustment parameters are programmable to take care of any minor mismatch of different sub assemblies & Balances are configured with these parameters.

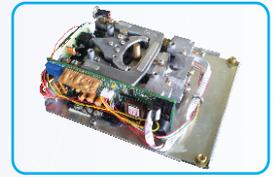

Advanced EMFC Transducer

Flameproof cum weather proof enclosure with glass window, push button, rotary switch in cast aluminium construction in compliance with requirements in respect of flameproofness and temperature class suitable for use in Zone 1 & Zone 2 and gas groups IIA & IIB atmospheres. Certified by Central Institute of Mining & Fuel Research, Govt. of India.

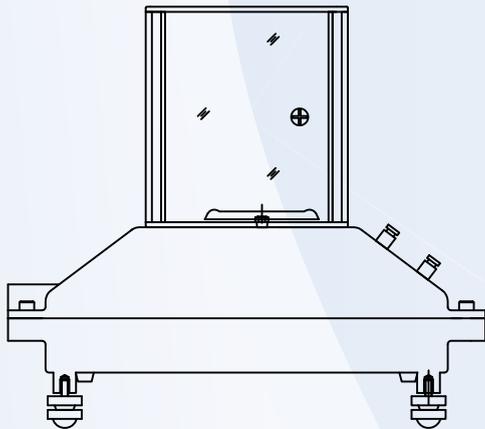


Range :

0.1 mg – 300g

1mg – 1000g

0.01g – 3kg

0.1g – 5kg

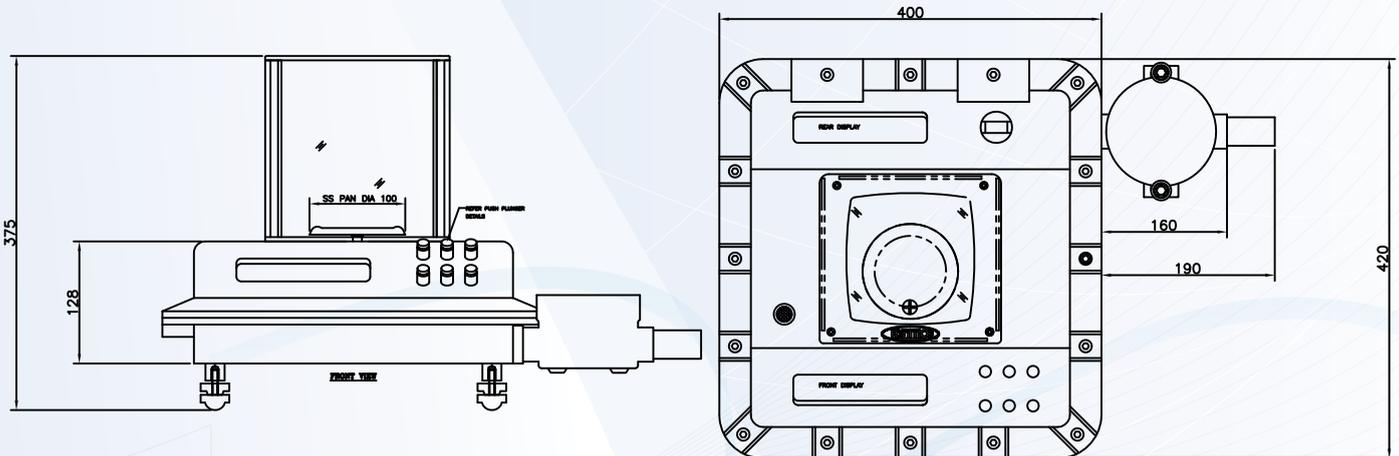

Specifications :

Taring Range : Up to full capacity

Response time : 3-5 seconds

Power Req : 90-250V, 50/60Hz, 10 Watt

Operating Temperature : 10 to 45 deg C

All Specifications and designs are subject to change for improvement.



CONTECH
Instruments Ltd.
AN ISO 9001, 14001 & OHSAS 18001 COMPANY

Where innovation never stops.....

301, Punit Indl. Premises, Turbhe Naka, Navi Mumbai - 400 705.

Tel. : +91 22 6139 3000 (100 lines)

Fax : +91 22 2761 8421

E-mail : info@contechindia.in / sales@contechindia.in

Website: www.contechindia.com



State Of The Art Manufacturing Plant